

myCobot Pro force-controlled gripper

1 Product Image



2 Specifications

name	myCobot Pro force-controlled gripper
Material	PC、PBT
size	156X106X61mm
Process Technology	Injection Molding
Gripping range	0-100 mm (default fingertip)

name	myCobot Pro force-controlled gripper
Repeatability accuracy	0.5 mm
Service life	300,000 opening and closing times
Drive Mode	Electric drive
Transmission mode	Gear + connecting rod
size	158x105x55mm
weight	340 g
Rated load	500g
Operating voltage	24V
Fixation method	Screw fixing
Environmental requirements	Normal temperature and pressure
Control interface	RS485/IO control/button control
Applicable devices	ER myCobot 320 series,ER Mercury series,ER myCobot Pro 600,ERmyCobot Pro 630,other general purpose robots

3 Working Principle


Driven by the motor, the finger surface of the manipulator makes linear reciprocating motion to achieve opening or closing. By setting the clamping torque, the impact on the workpiece is minimized, the positioning point is controllable, and the clamping is controllable.

4 Usage scenarios

Experimental operation: In scientific research experiments, complete the grabbing and moving of test tubes and utensils to ensure the safety and accuracy of the experiment. Educational demonstration: As a teaching tool, it helps students understand the principle of robot grabbing and cultivate practical ability. Material handling: In simulated production lines or warehouses, materials of various specifications are moved to improve work efficiency.


5. Installation

Install the gripper connector to the end-of-arm flange using screws and washers

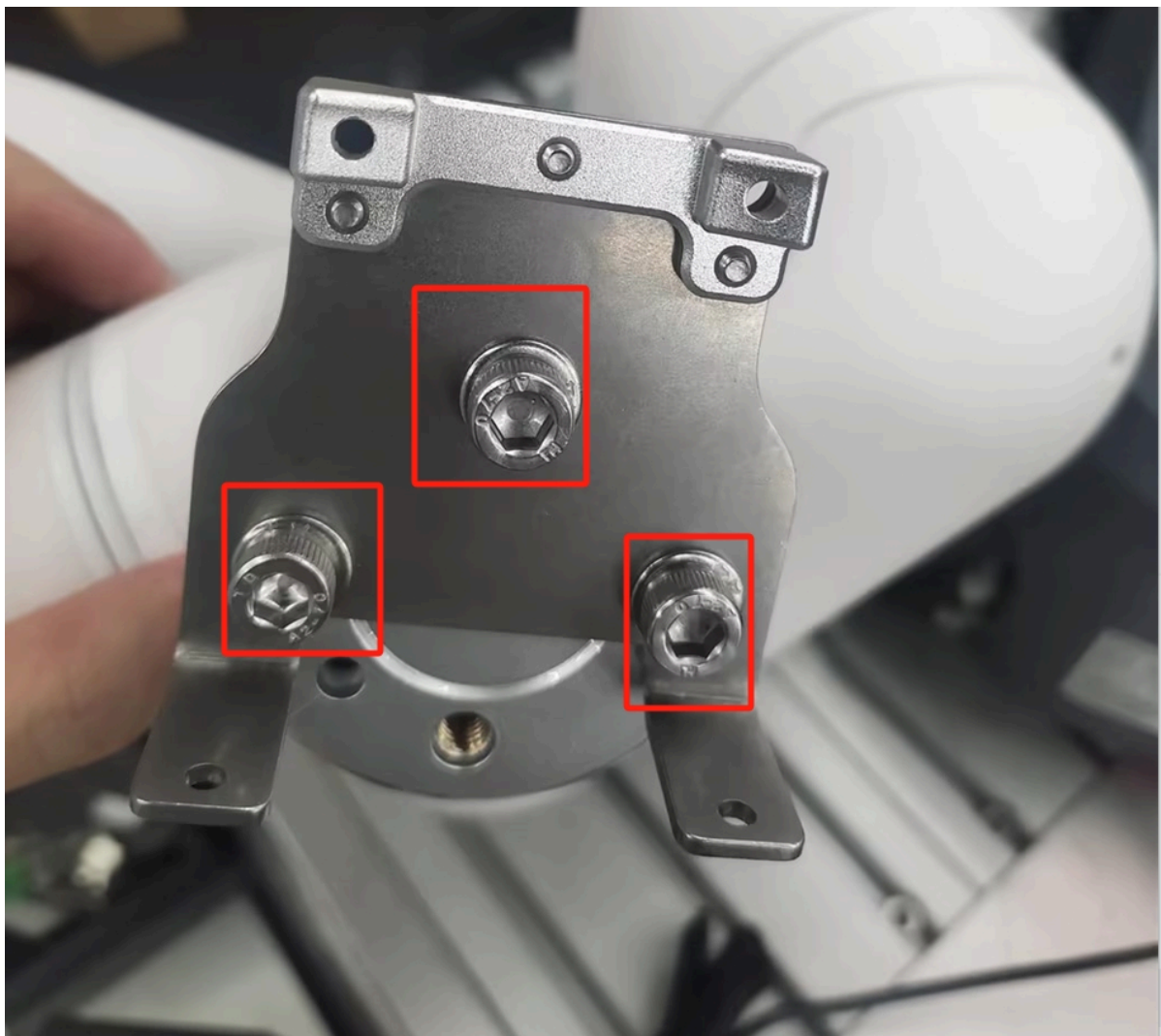

Then use screws to install the clamp on the connecting piece.



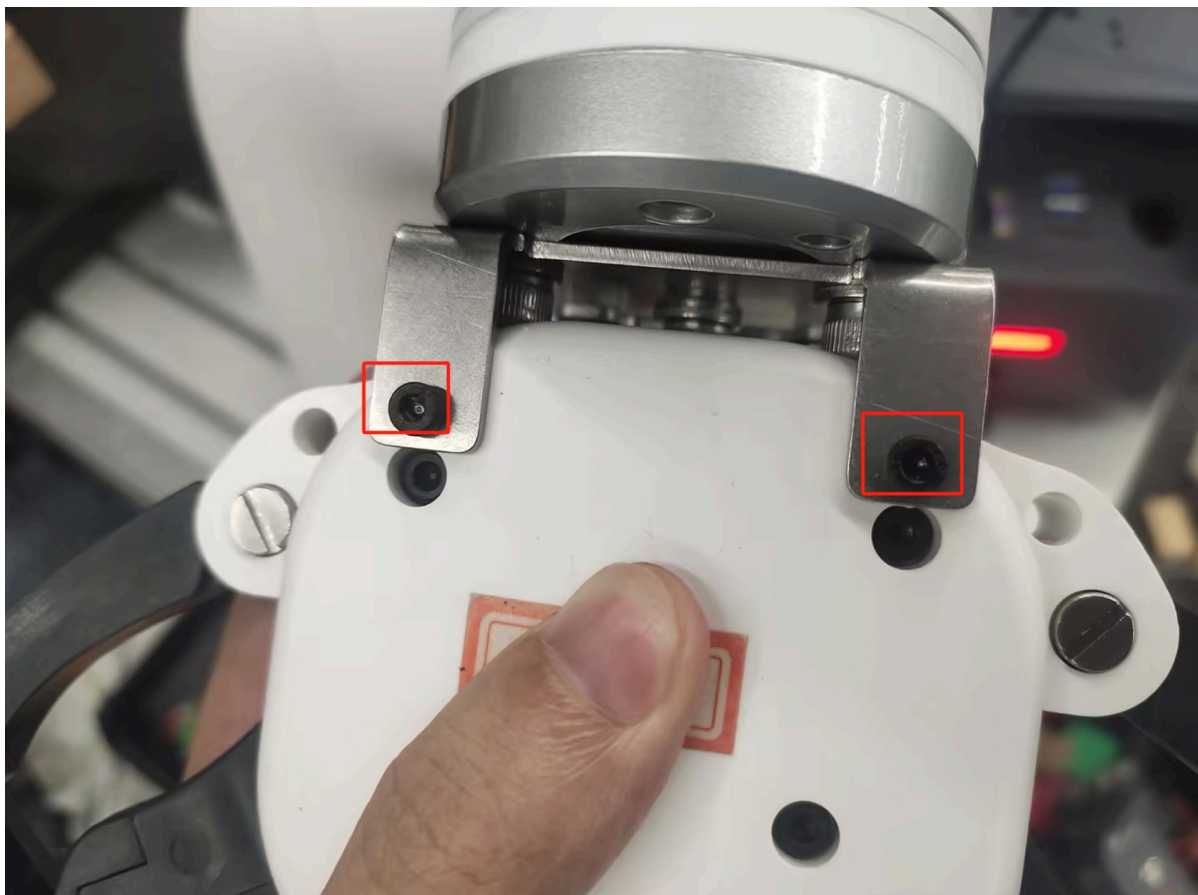

Finally, use M8 aviation wire to connect the gripper and the robotic arm.

